

PremiCron®

Non-absorbable sutures


DESCRIPTION

PremiCron® is a non-absorbable braided and coated synthetic multifilament composed of polyester fibres. PremiCron® sutures are coated uniformly with silicone to increase surface lubricity, thereby enhancing handling characteristics.

Available in single packs (1 or 2 sutures) or multipacks (4 or 8 sutures) either green or white.

INDICATION

- Cardiac surgery (valve surgery)
- General surgery

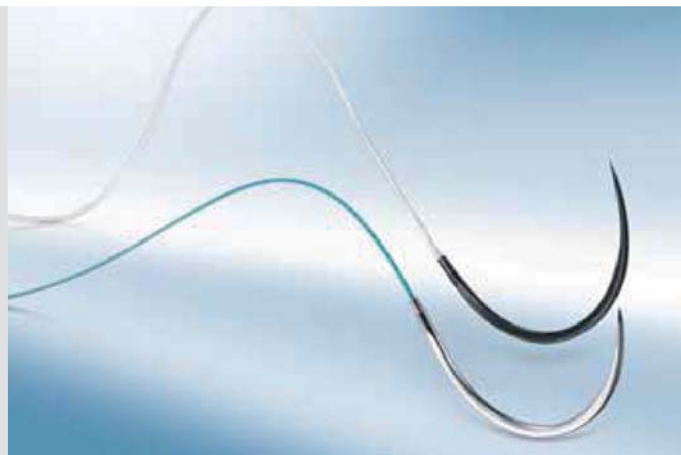

PremiCron®

ADVANTAGES

- PTFE oval pledget design improving adaptability to underlying anatomical structures
- Excellent knot run-down and knot security
- Reduced tissue drag and sawing
- 300 series stainless steel needles for superior breaking resistance "Easyslide"
- New generation of needles for cardiovascular indication: CV pass and CV pass easyblack

PremiCron®

Product overview

Structure

Braided



Colour

Green or white

Chemical composition

Polyester

Coating

Silicone

Origin

Synthetic

Sizes

USP 6/0 (metric 0.7) to USP 5 (metric 7)

Type of absorption

Non-absorbable

Sterilisation

Ethylene oxide / Gamma-irradiation